


Increase liquid handling quality with easy and reliable performance verifications.




Artel MVS®

Without knowing the exact volume transferred at each step of an assay, the concentration of species in solution is unknown, potentially leading to unreliable results, delays, and cost overruns. Utilization of the Artel MVS for performance verification can prevent problems, bolster laboratory productivity, and provide confidence in the integrity of data.


Three simple steps to volume verification with automatic reporting.

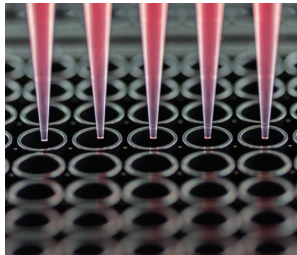

1. Dispense the Sample Solution and Diluent into a Verification Plate and mix.

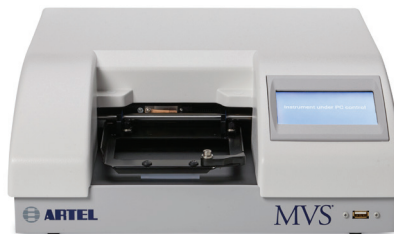

2. Measure the absorbance values with the Plate Reader.

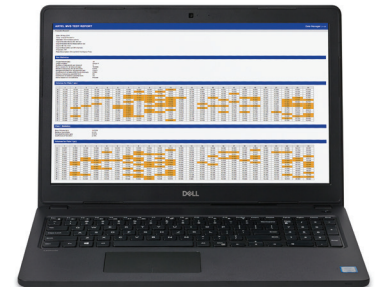

3. Software displays results that indicate the performance of the liquid handler.

Ideal for:

- Verifying the accuracy and precision of automated liquid handlers and multichannel pipettes
- Comparing instrument performance regardless of make, model or location
- Obtaining highly accurate results traceable to national and international standards
- Standardizing results among multichannel pipette operators

Verify volume transfer precision and accuracy simply and traceably in your lab, and simplify regulatory compliance with the easy-to-use MVS.

System Capabilities:

- Measures performance over a wide volume range of 0.0001 to 350 μ L for Automated Liquid Handlers and Multichannel Pipettes
- Measures aqueous and non-aqueous volume transfers, including dimethyl sulfoxide (DMSO)
- Obtains volume statistics tip-by-tip and well-by-well
- Eliminates the need for rigorous environmental controls
- Simplifies verification and calibration—easy-to-use for technicians of any skill level

Key Applications:

- Facilitating of scale-up and method transfer
- Developing, troubleshooting, and validating assays, methods, and SOPs
- Ensuring reproducibility
- Site acceptance testing and establishing baseline performance of new equipment
- Verifying performance before and after maintenance
- Assessing performance of other essential liquid handling equipment such as plate shakers, plate washers, and bulk dispensers
- Assessment and training of multichannel pipette operators



“The MVS is easy to use, facilitating a frequent performance monitoring schedule. has made a big difference in our laboratories, both from an efficiency point of view, as well as from a confidence standpoint.”



artel-usa.com/MVS